


Collection Protocol


Deposit saliva into the **Collection Tube** up to the indicated fill line.

If you notice excess bubbles, eliminate them by gently tapping the tube on a hard surface.

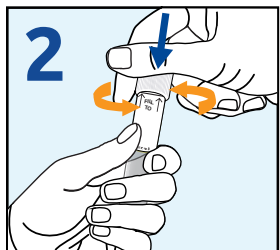


Remove the **Collection Funnel**. **Cap** the Collection Tube by **pushing down** and **twisting** until secured, and the **Yellow Saliva DNA Preservative** will be released.

Discard or recycle the Collection Funnel.

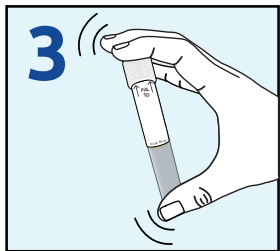

Shake well for 10 seconds to mix the saliva with the **Yellow Saliva DNA Preservative**.

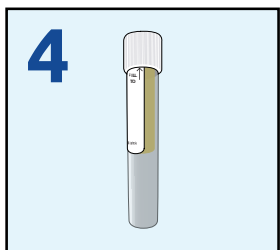

The **Collection Tube**, containing the saliva sample is ready for transport.

Intended Use:

Norgen's EZCollect Saliva DNA Device is designed for 1) simple and non-invasive home-based saliva collection, and 2) preservation of DNA in saliva samples at ambient temperature. Intended for Research Use Only.

Storage:

This device should be stored at room temperature for up to 2 years without any reduction in performance prior to Saliva Sample Collection. After the Saliva Sample Collection, the saliva DNA in preserved samples is stable for up to 6 years at room temperature (15-30°C).

Before starting the sample collection:

Do not eat, drink, smoke, chew gum, brush or floss for 30 minutes before collection. Rub your cheeks and wiggle your tongue to stimulate the secretion of saliva. It may take up to 7 minutes to collect the sample. Do not exceed 30 minutes for sample collection.

Warnings:

If any of the Saliva DNA Preservative comes into contact with skin or eyes, wash thoroughly with water. Do not ingest the Saliva DNA Preservative. The collection device contains a choking hazard (Collection Tube Cap). Keep away from children. The Saliva DNA Preservative and preserved saliva samples will form reactive compounds when combined with bleach, thus care must be taken to properly dispose of any spilled preservative or remaining samples.

Tech support: techsupport@norgenbiotech.com
Orders: orders@norgenbiotech.com

Made in CANADA

Rev. 1



Scan for instructional video or visit
norgenbiotech.com/ezcollect/salivadna

Patent Pending