

Solutions for small volumes.



The OD600 DiluPhotometer™



The OD600 DiluPhotometer™

Virtual Dilution Built-in

Small, portable spectrophotometer for quantifications at 590 to 610 nm
Configured for bacterial growth rates and colorimetric protein assays

MAINTENANCE FREE SETUP

No moving parts, no recalibration necessary

BATTERY OPERATED

Unlimited mobility, ideal for usage in incubation cabinets and under anaerobic conditions

DATA OUTPUT VIA RS232 or USB INTERFACE

Submit data on a computer

BUILT-IN MEMORY FOR DATA STORAGE

Storage of up to 99 measurements

GLP/GMP READY WITH DILUCELL™


Small volume requirement 200 μ l or 100 μ l
Automatic sample dilution of factor 10 and 20

FAST AND EASY HANDLING


Measurement in less than 3 seconds




OD600 Application

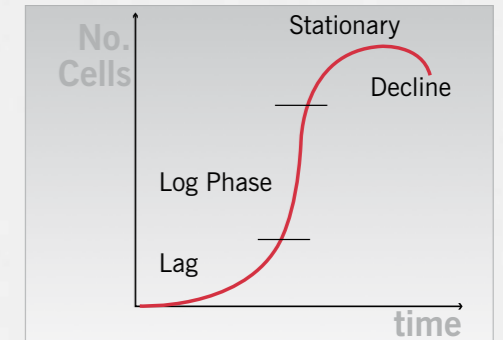

Protein Application

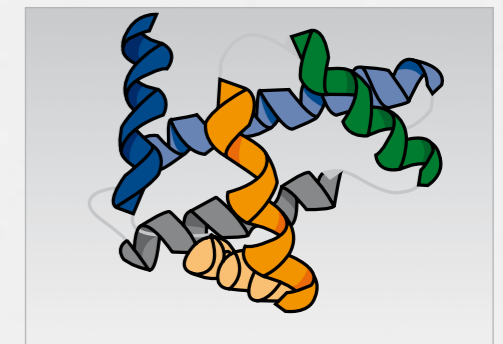

Virtual Dilution with DiluCell™

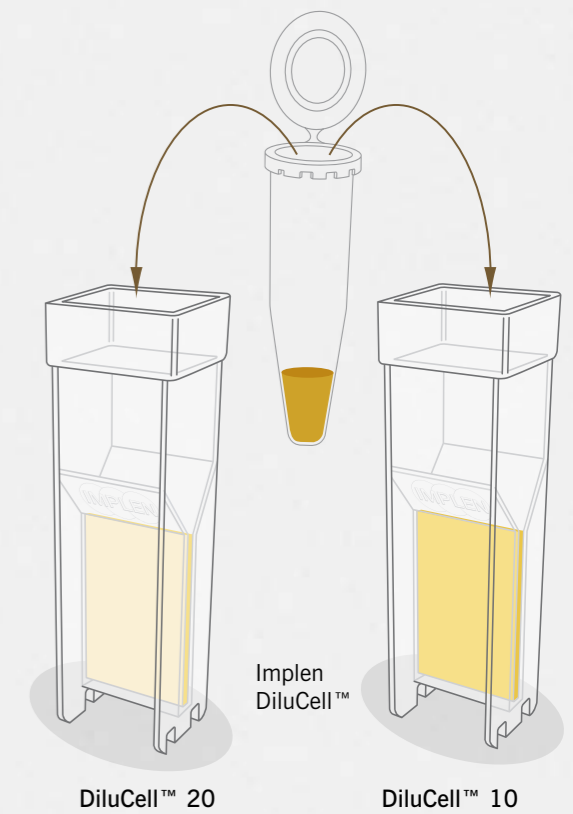

Technical Specifications

OD600 DiluPhotometer™

Wavelength	600 nm
Bandwidth	40 nm
Range	Optical Density – 0.3 A to 1.99 A
Accuracy	< ±0.05 A at 1 A using Neutral Density Filters
Repeatability	±0.02 A at 1 A
Controls	6 button – power, blank, sample, memory, reset, recall/print
Optical height	8.5 mm
Cuvette holder	Fixed with drain hole. Accepts DiluCell™ 10, DiluCell™ 20, 10 mm pathlength semi micro and macro cuvettes or 14 – 16 mm round tubes.
Output	RS232, USB
Memory	99 readings
Display	Custom LCD
Power requirements	External power adaptor (110 to 220 V, 50/60 Hz, 20 VA) or internal rechargeable NiMH battery
Dimensions	180 x 150 x 60 mm
Weight	0.6 kg

DiluCell™

Dimension	12 mm x 12 mm x 45 mm
Pathlength	0.5 mm (DC 20) or 1.0 mm (DC 10) (identifyable by special marking)
Automatic dilution	factor 20 (DC 20) or factor 10 (DC 10)
Volume	DC 20: 100 µl max. sample volume DC 10: 200 µl max. sample volume
UV/VIS transparent	starting from 340 nm
Availability	Approved for usage with the OD600 DiluPhotometer™

Features and specifications are subject to change without notice.



Europe and Asia

Implen GmbH
Schatzbogen 52
81829 München
Germany

Phone +49 89 72637180
Fax +49 89 726371854
Email: info@implen.de
www.implen.de

North and South America

Implen, Inc.
Suite 521 – 23
5655 Lindero Canyon Road
Westlake Village
CA 91362
USA

Phone +1 818 706-3500
Fax +1 818 449-0416
Email: info@implen.com
www.implen.com

UK Office

Implen UK, Ltd.
Insight House
Riverside Business Park
Stoney Common Road
Stansted Mountfitchet
Essex, CM24 8PL
Great Britain

Phone +44 1487 843722
Fax +44 1480 260111
Email: info@implen.co.uk
www.implen.co.uk

Switzerland Office

Implen Schweiz
Schwaentenmos 15
CH-8126 Zumikon
Switzerland

Phone +41 44 9180896
Fax +41 44 9180812
Email: info@implen.ch
www.implen.ch